


SparMED Aps
Ryttermarken 2
3520 Farum
Denmark

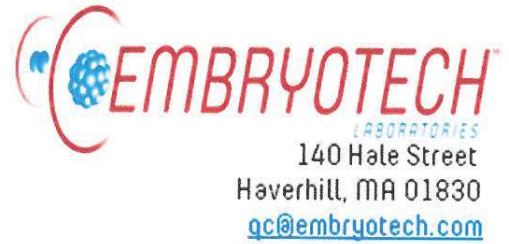

ELI Accession Number: SPAR-6493-0317

Date of completion: 03-26-2017

Lot number: 07730

Reference number: OOPW-SW02

Description of test article(s): Oosafe® 6-Well Dish with Straw Holder

Assay system requested by customer: 0.25mL of culture medium was placed in each well of the test article and overlaid with oil. 3-4 one-cell mouse embryos were added to each of the six wells and cultured for 96-hours.

Control assay method and results: 15 one cell (B₆C₃F₁ X B₆D₂F₁) embryos were cultured in a 1mL drop of culture medium:

15 / 15 (100 %)

1-cell to 2-cell within 24 hr

15 / 15 (100 %)

1-cell to expanded blastocyst within 96 hr

For a valid assay, Embryotech™ requires at least 70% of one cell stage control embryos to develop to expanded blastocyst within 96-hours.

Test assay method and results: 21 one cell (B₆C₃F₁ X B₆D₂F₁) embryos were cultured in the test article using culture medium:

21 / 21 (100 %)

1-cell to 2-cell within 24 hr

21 / 21 (100 %)

1-cell to expanded blastocyst within 96 hr

Summary of observations: All test and control embryos were selected randomly from a common pool of freshly collected embryos and were cultured in the same incubator at 37°C and 5.0% CO₂. 100 percent of the control embryos developed to the expanded blastocyst stage within 96-hours. 100 percent of the embryos cultured in the test article developed to the expanded blastocyst stage within 96-hours.

Signature
Study Director

03-27-2017
Date

Signature
Quality Reviewer

03-27-2017
Date