




SparMED Aps
Ryttermarken 2
3520 Farum
Denmark

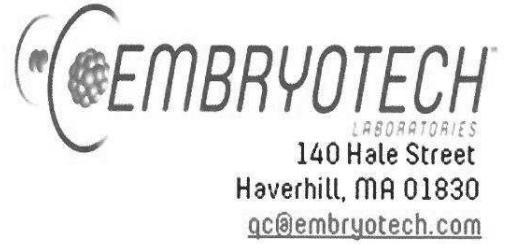

140 Hale Street
Haverhill, MA 01830
qc@embryotech.com

ELI Accession Number: SPAR-8274-0218

Date of completion: 02-11-2018

Lot number: 07848

Reference number: OOPW-FW03

Description of test article(s): Oosafe® 4 Well Dish - Treated Surface

Assay system requested by customer: 0.5mL of culture medium was placed in the each well of the test article and overlaid with oil. 5-6 one cell mouse embryos were added to each well of the test article and cultured for 96-hours.

Control assay method and results: 15 one cell (B6C3F1 X B6D2F1) embryos were cultured in triplicate micro drops of culture medium:

15 / 15 (100 %)

1-cell to 2-cell within 24 hr

15 / 15 (100 %)

1-cell to expanded blastocyst within 96 hr

For a valid assay, Embryotech™ requires at least 70% of one cell stage control embryos to develop to expanded blastocyst within 96-hours.

Test assay method and results: 21 one cell (B6C3F1 X B6D2F1) embryos were cultured in each well containing culture medium overlaid with oil in the test article:

21 / 21 (100 %)

1-cell to 2-cell within 24 hr

21 / 21 (100 %)

1-cell to expanded blastocyst within 96 hr

Summary of observations: All test and control embryos were selected randomly from a common pool of freshly collected embryos and were cultured in the same incubator at 37°C and 5.0% CO₂. 100 percent of the control embryos developed to the expanded blastocyst stage within 96-hours. 100 percent of the embryos cultured in the test article developed to the expanded blastocyst stage within 96-hours.

Signature
Study Director

02-12-2018
Date

Signature
Quality Reviewer

02-12-2018
Date