

# Certificate of Analysis

**Product name:** Andrology Tube, 5 ml  
**Product type:** Oosafe® Plasticware IVD  
**Product code:** OOPW-AT10  
**Batch/LOT No.:** 09034  
**Product expiry date:** 2029-10

TESTING METHOD	ACCEPTANCE CRITERIA	RESULTS	CONCLUSION	DOCUMENT NO.
Mouse Embryo Assay Test (MEA)	Over 80% embryo development to the expanded blastocyst stage within 96 hours	Passed	Non-embryotoxic	SPAR-5130-1024
Human Sperm Survival Assay (HSSA)	≥75% sperm motility after 24 hours proven	Passed	Stable human sperm motility	S8289-1024SPAR
Limulus Amoebocyte Lysate (LAL)	<0.03 EU/device	Passed	Do not contain endotoxin	E4290-1024SPAR

Quality Management system according to ISO 13485:2016

**Sterilized by irradiation according to ISO 11137-1**

**Sterilization method:** E-Beam

**Sterility assurance level:** 10<sup>-6</sup>


**GOosafe with SparMED!**

Manuela Jagieniak

Quality Control Department

SparMED ApS

Issue date: 29.10.2024

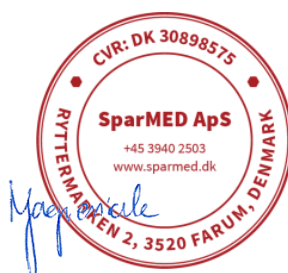

**Sparmed ApS**

Ryttermarken 2, 3520 Farum, DK ● +45 39 40 25 03 ● info@sparmed.dk ● www.sparmed.dk  
 Corp. VAT no.: DK 30 89 85 75