

PREFACE

For our last 2021 Gynemedia we wanted to talk about our GYNEMED products: GM501 Mineral Oil, paraffin oil ready to use, and our Oocyte aspiration needles.

Also share will you are participation to the MEFS conference in Dubai, our first international event since the beginning of the pandemic!

And finally to inform you about our close days for Christmas.



GYNEMED GM501 Mineral Oil

GM501 Mineral Oil is a ready-to-use oil for covering the medium during IVF/ICSI treatment.

Key features:

- Stabilize the pH and temperature during IVF, ICSI and related artificial reproductive techniques
- CE marked class IIa (0482)
- MEA test (one-cell assay, 120 h exposure)
- Filled in glass flasks to avoid peroxide built-up
- Available in 100mL and 500mL
- Shelf Life 18 month from time of manufacture

Composition:


- Paraffin oil
- Density 0.83 - 0.86 g/ml
- Viscosity < 30 cP at 30°C
- Pre-washed twice with ultra-pure water



Overlay the culture medium with GM501 Mineral Oil until it is completely covered. We recommend preparing the dishes with oil overlay the day before use. If using bicarbonate buffered media incubate the dishes in a CO₂ incubator at 37°C to ensure warming and saturation with gas. Each laboratory needs to establish and validate its own procedures for medium preparation and use.



Oocyte Aspiration Needle

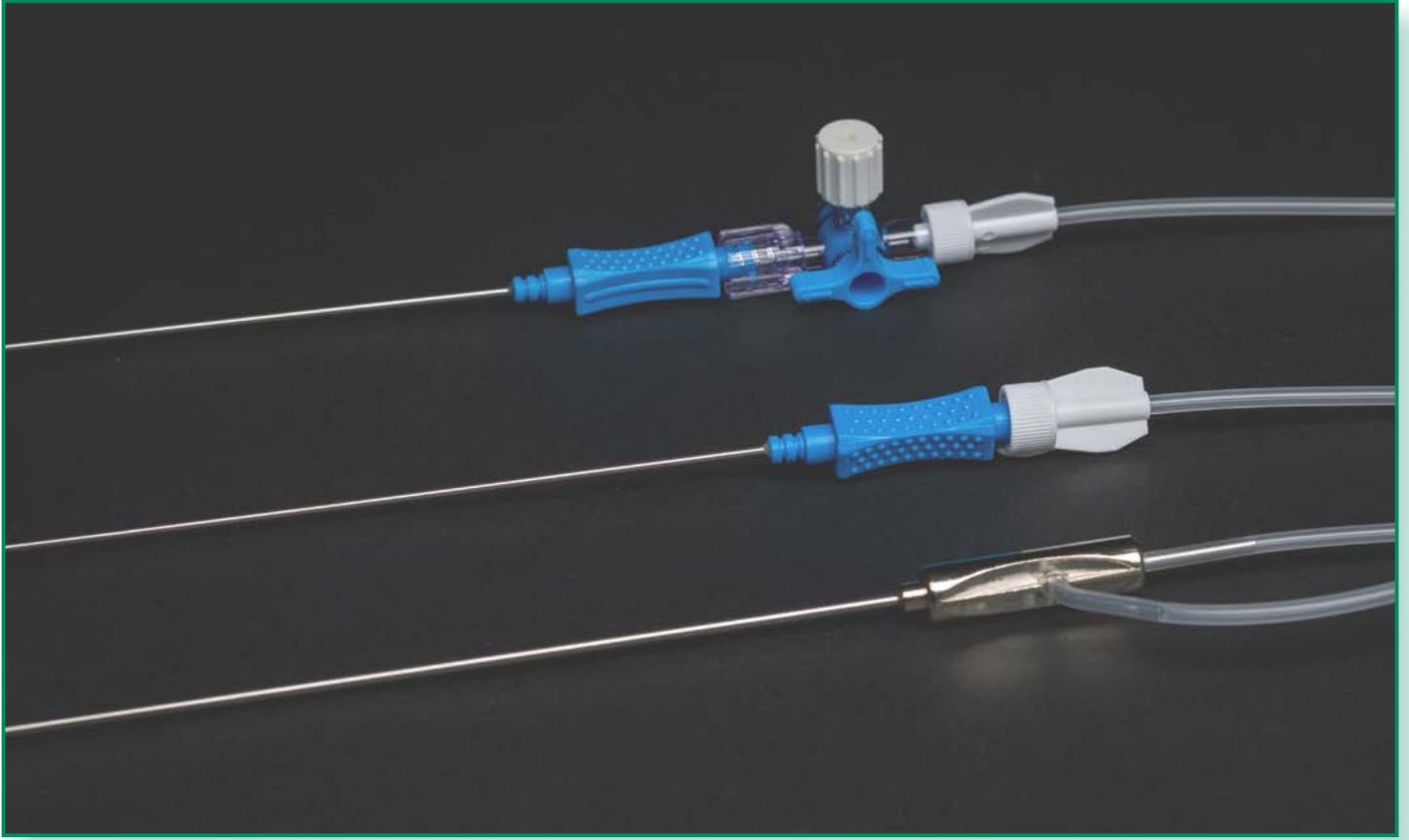

The oocyte aspiration needles are available as single or double lumen sets.

The single lumen needles are available with or without flushing

(3-way-valve). The double lumen needles can be used for concurrent aspiration and flushing (stopcock directly connected to handle). The silicone stopper is suitable for 14 ml-round bottom collection tubes.

The oocyte aspiration needles have comfortable handles with position marking.

Variable needle and tubing configurations are available (length, thickness etc.).

Key features:

- **Sterile** - Sterilized by gamma irradiation and individually packed.
- **Echogenic markings** - Enable excellent visualization under ultrasound view.
- **Customized production** - Available upon request.
- **CE marked** - To meet international standards (Class IIa).
- **Shelf life** - Sterility and functionality guaranteed for 3 years from time of manufacture.
- **MEA-tested** - A certificate of analysis (COA) is available for each lot number upon request from our website.
- **Packed by 10 units**

Available at Gynemed!

MEFS - Middle East Fertility Society

In December you can meet us at the MEFS in Dubai, taking place at the InterContinental Festival City.

After two years we are back in this region to meet with our customers!

At our booth we will be presenting our GYNEMED product line but also Hamilton Thorne equipment and Planer products.

Come to visit us at the booth number 6

Our international sales manager will be on site, do not hesitate to contact her

laurence.rivain@gynemed.de

2nd-4th of December 2021

www.mefs.org



CHRISTMAS CLOSURE



Gynemed will be closed from the 24th until the 31st of December, included.
Don't forget to anticipate your order during this period!
Last day of shipment is 21st of December.

LEGAL NOTE

Publisher: GYNEMED GmbH & Co. KG
Telephone: +49 4363/90329-0 Fax: +49 4363/90329-19 Email: info@gynemed.de
Editor: Dr Julia Heinzmann (V.i.S.d.P.) Layout: Julia Biegemann